


# INTEOS® - Radius Plate palmar

- watershed line design - distal eXtreme
- watershed line design - standard

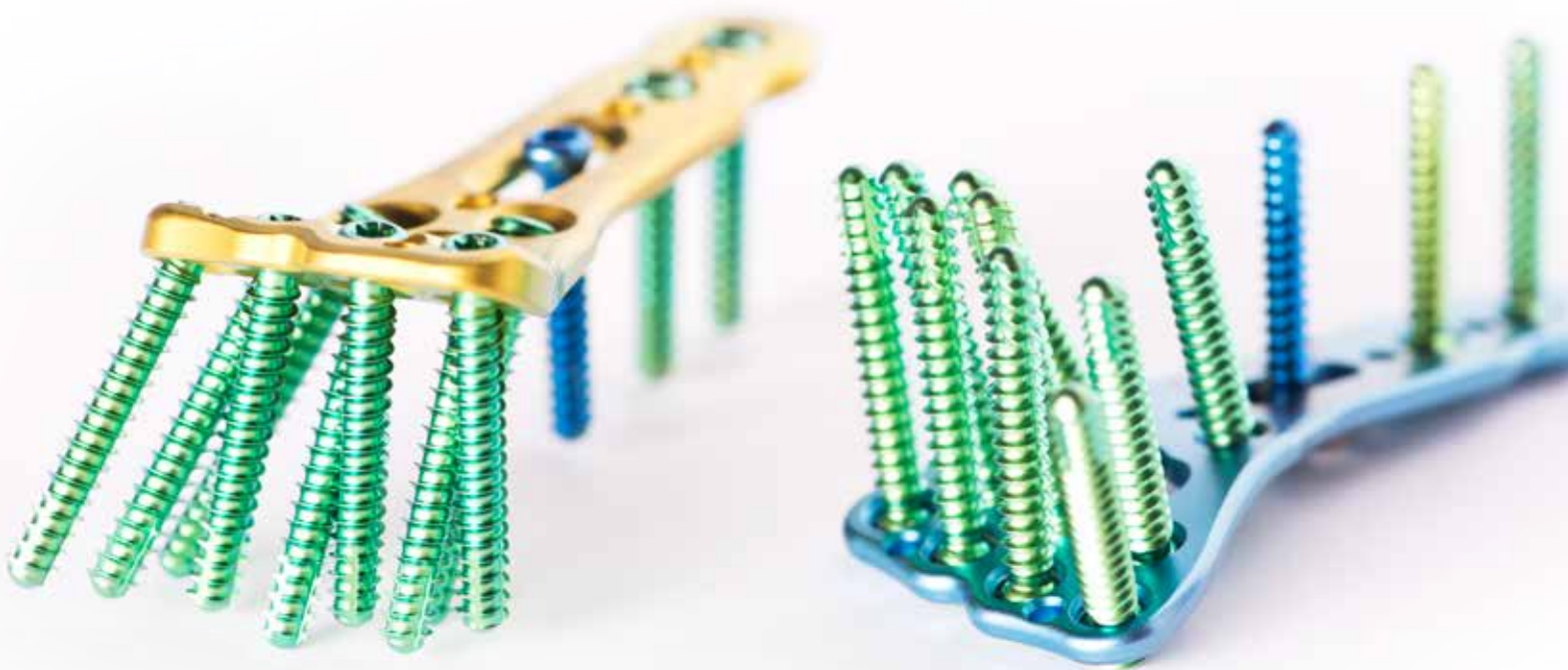

# Advantages of the INTEOS® Radius Plate wsl & pal

The INTEOS® Radius Plate series M3 of Hofer-medical are characterized by its easy implantation technique as a result of pre-angled holes. **Fixed angle & multidirectional** implantation techniques come as standard.

INTEOS® - Technology




## Plate features

- Minimal risk of perforation for placement close to the joint as a result of pre-angled, anatomically optimized holes
- Increased dorsal stability in the osteoporotic bone thanks to a third (supporting) row of screws
- Anatomically optimized for the ulnar key fragment
- Consistent watershed line design
- Minimal screw head protrusion and overlap
- Various possibilities of reduction via reduction instrument through K-wire holes and windows in the plate
- Direct fixing of the styloid to the plate

## Screw features

- Rounded screw head
- Non-invasive screw tip
- „Stick fit“
- ForceDrive 8

Fixed angle and multidirectional functionality thanks to the INTEOS® Small Fragment System standard

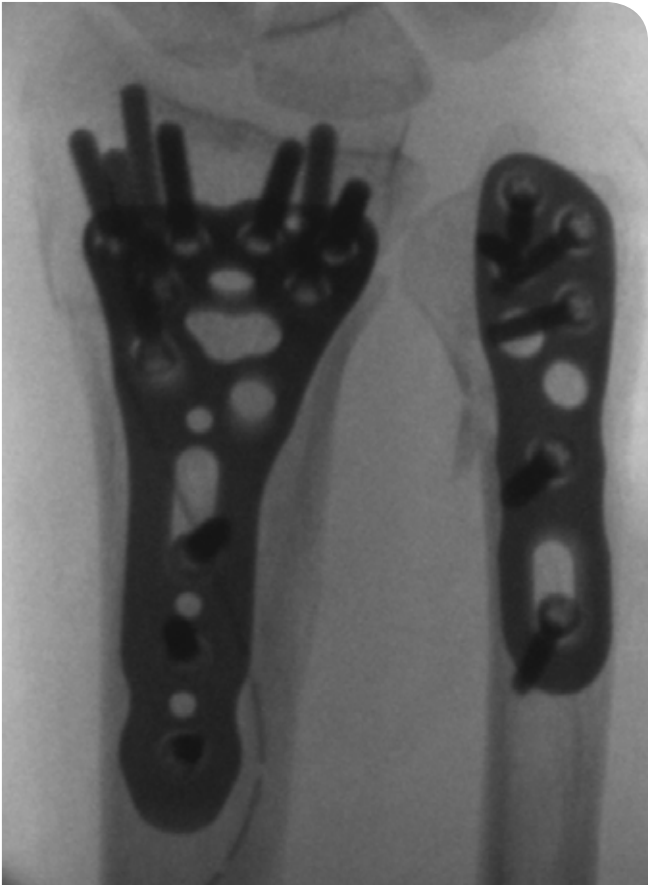


Watershed line design  
Distal eXtreme

Watershed line design  
Standard



Watershed line design as well as perfect anatomical fit of the system

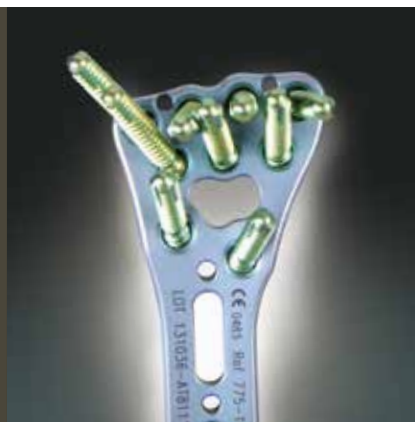

Dorsal stability in the osteoporotic bone owing to a third row of screws



Innovative screw placement to support the joint



**AUSTRIA**

HOFER GMBH & CO KG  
Jahnstrasse 10-12  
A 8280 Fürstenfeld  
TEL: +43(0)3382/53388  
FAX: +43(0)3382/53093  
office@hofer-medical.com  
www.hofer-medical.com  
FN21826y  
UID ATU30764704



**GERMANY**

HOFER-MEDICAL GMBH  
Ruhrtalstraße 52-60  
D 45239 Essen  
TEL: 0800 -700 80 22 (innerhalb BRD)  
FAX: +49 201 56599 1069  
office@hofer-medical.de  
www.hofer-medical.com  
HRB24645  
UID DE290596924



**ITALY**

HOFER MEDICAL ITALIA S.R.L.  
Via della Giuliana, 66  
I 00195 Rome

TEL: +39 393 8950481  
FAX: +39 06 81151148  
office@hofer-medical.it  
www.hofer-medical.it  
RM-1454597  
UID IT13534721009



**POLAND**

HOFER GMBH & CO KG SPÓŁKA  
KOMANDYTOWA ODDZIAŁ W POLSCE  
Radkowska 9  
PL 57-402 Nowa Ruda

TEL.Kom.:+48 603224112  
FAX: +48 748712273  
office@hofer-medical.pl  
www.hofer-medical.com  
REGON 141527293  
UID PL1070011771

## creating partnership

For us, a partnership is the result of a long and happy relationship in all areas of our work. Reliable bone healing for patients, optimal handling of the systems for surgeons and surgical personnel, as well as the simple preparation of the instruments, constantly strengthen this partnership.

For us, "creating" means more than just finding solutions in the form of innovative products.

It is the result of high standards, constant development, innovative products and excellent service for patients, surgeons and surgical personnel around the clock.

Please do not hesitate to contact us if you have questions about our company, our employees or our production methods.

Your HOFER medical solutions Team

vCard

